F

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E

## PACKAGE INFORMATION LEAFLET

**Aktimed GRID TAPE CLASSIC** is a lattice-based tape based on polyurethane with an acrylic adhesive.

### APPLICATIONS

- Pain Point, Trigger Point and Acupuncture Point Taping\*)

### INDICATIONS FOR USE

- The product is intended for single use.
- Application is recommended by a trained person such as a doctor, naturopath or physiotherapist.
- Before applying the mesh plaster degrease the skin.
- Glue loses its adhesive power by heat.
- The tape should be applied 20-30 minutes before physical activity.
- It may occasionally cause skin irritation. In this case, immediately remove the tape.

### SAFETY INFORMATION

- In case of allergy to any of the ingredients do not use the tape or first test on a small area of skin!
- Do not adhere to open wounds or skin injuries!
- In the case of skin reactions such as itching or redness remove the tape immediately!
- Pregnant women may only apply the tape after consultation with the attending physician!
- Keep out of reach of children!

\*) Effects have not been scientifically proven

E

### APPLICATION TIME

The tape lasts up to 6 days on the skin.

### INGREDIENTS

80% polyester, 15% polyurethane, 5% acrylate.

### TECHNICAL DATA

One pack contains 20 leaves each with different numbers of tapes depending on their size:

- **SMALL:** 2.1 cm x 2.7 cm - one sheet contains 9 tapes
- **MEDIUM:** 2.8 cm x 3.6cm - one sheet contains 6 tapes
- **LARGE:** 4.4 cm x 5.2cm - one sheet contains 2 tapes

Store in a dry and cool place away from sunlight.  
Store between 8° C and 25° C.


CE: conforms to EU Directive 93/42/EEC for medical products.

Aktimed GmbH, Theaterstr. 18, 69117 Heidelberg, Germany  
www.aktimed-tape.de

**Aktimed**  
GRID TAPE CLASSIC

**GRID TAPE**  
CLASSIC

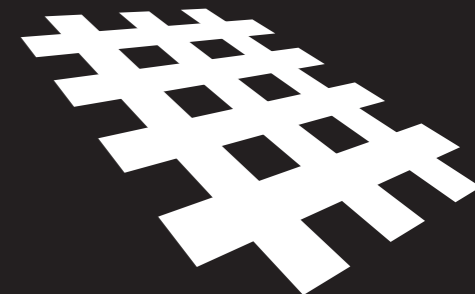D E F I